


4.2 Perform the installation work at the customer's.

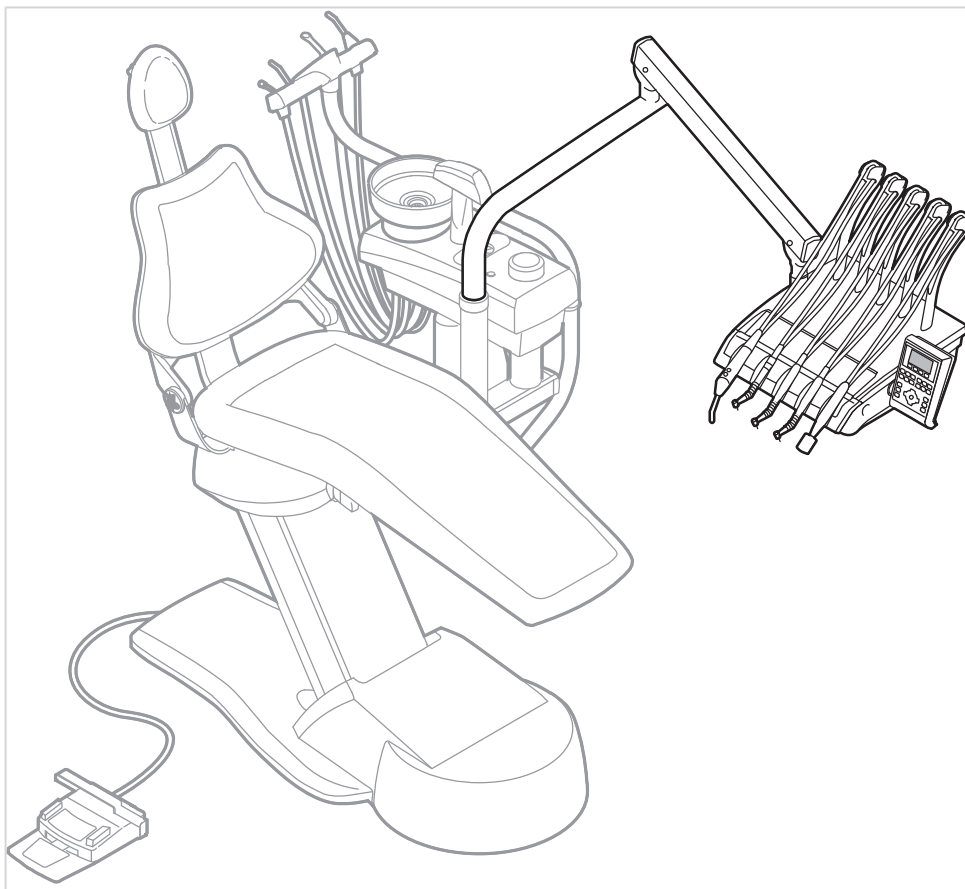

ESTETICA E30 S

Documents:

- Set-up plan
- Drilling template paper (**Mat. no. 1.009.2781**)
- Technical specifications and requirements


The connection diagram for the customer's installation is drawn on the drilling template.

All information on surface-mounted connection, flush-mounted connection and the floor attachment points are shown on the drilling template.

The drilling template has been drawn on a scale of 1:1 and can be used as a fitting aid.

- ▶ Identify the setup site according to the setup plan.
- ▶ Connect the power, water, etc.

4.2.1 Flush-mounted connection

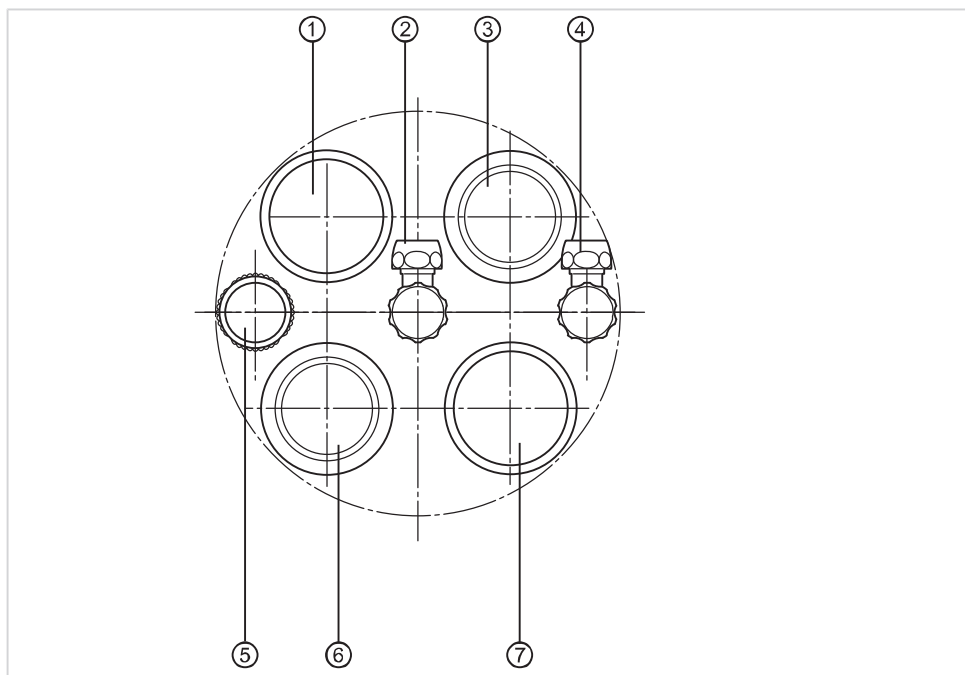

Flush-mounted connection diagram

- | | |
|-----------------------------------|----------------------------------|
| ① Connecting tube for multimedia | ⑤ Electric lines |
| ② Internal R 1/2 water connection | ⑥ NW40 drain connection |
| ③ NW 40 Suction connection | ⑦ Connecting tube for multimedia |
| ④ Internal R 1/2 air connection | |

Electrical lines:

- Electrical lead:
 - 3x1.5 mm² (customer fuse protection 10 A)
 - 3 x 2.5 mm² (customer fuse protection 16 A)
- Control line for central suction
- Protective conductor 1 x 4 mm²
- Bell line max. 24 V

Height of the air and water connections above the floor: min. 40 mm to max. 75 mm


4.2.2 Surface-mounted connection



Note

If the under floor connection is not located under the standard covering hood, an additional hood (**Mat. no. 1.009.5990**) can be ordered. The additional hood contains all the parts that are required for the assembly:

- Additional hood and holder
- Extension hoses for water and air connection
- Installation material

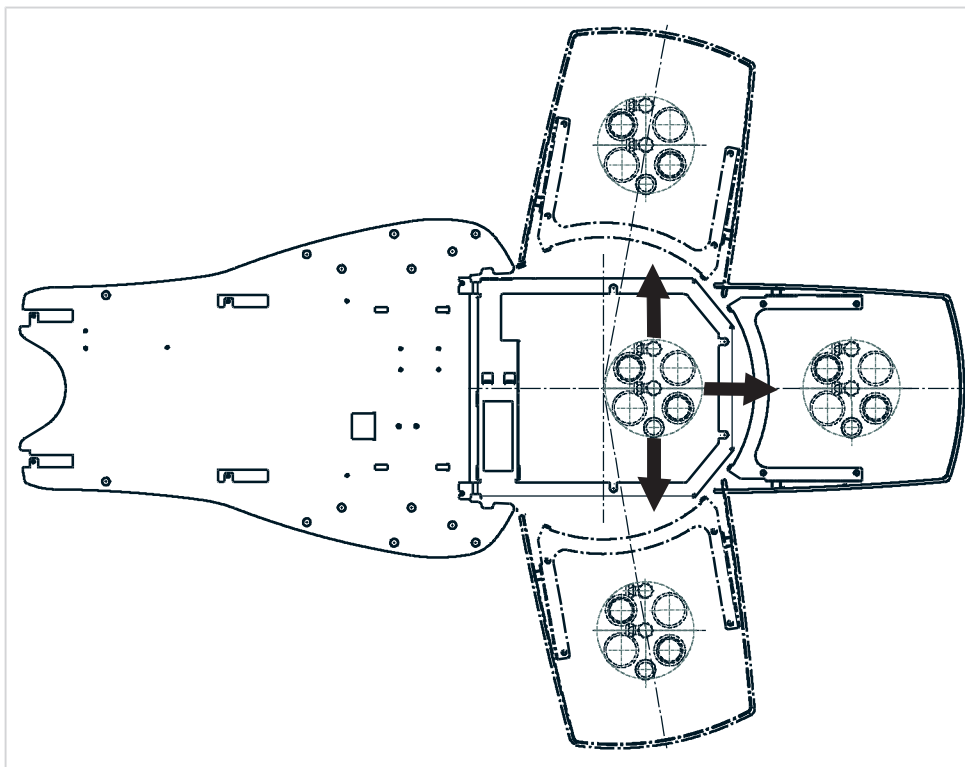

The additional hood can be installed from the central axis of the chair, in both directions and with a rotation of up to 80°.

- The detailed dimensions can be taken from the installation plan included in the delivery.

6 Mounting the kits




Note

Affixing and installing kits, please comply with the following:

- The kits may only be affixed and installed when the unit is shut off.
- For new installations, the kits must be mounted and installed before the device is commissioned.

7 Setup plans

7 Setup plans

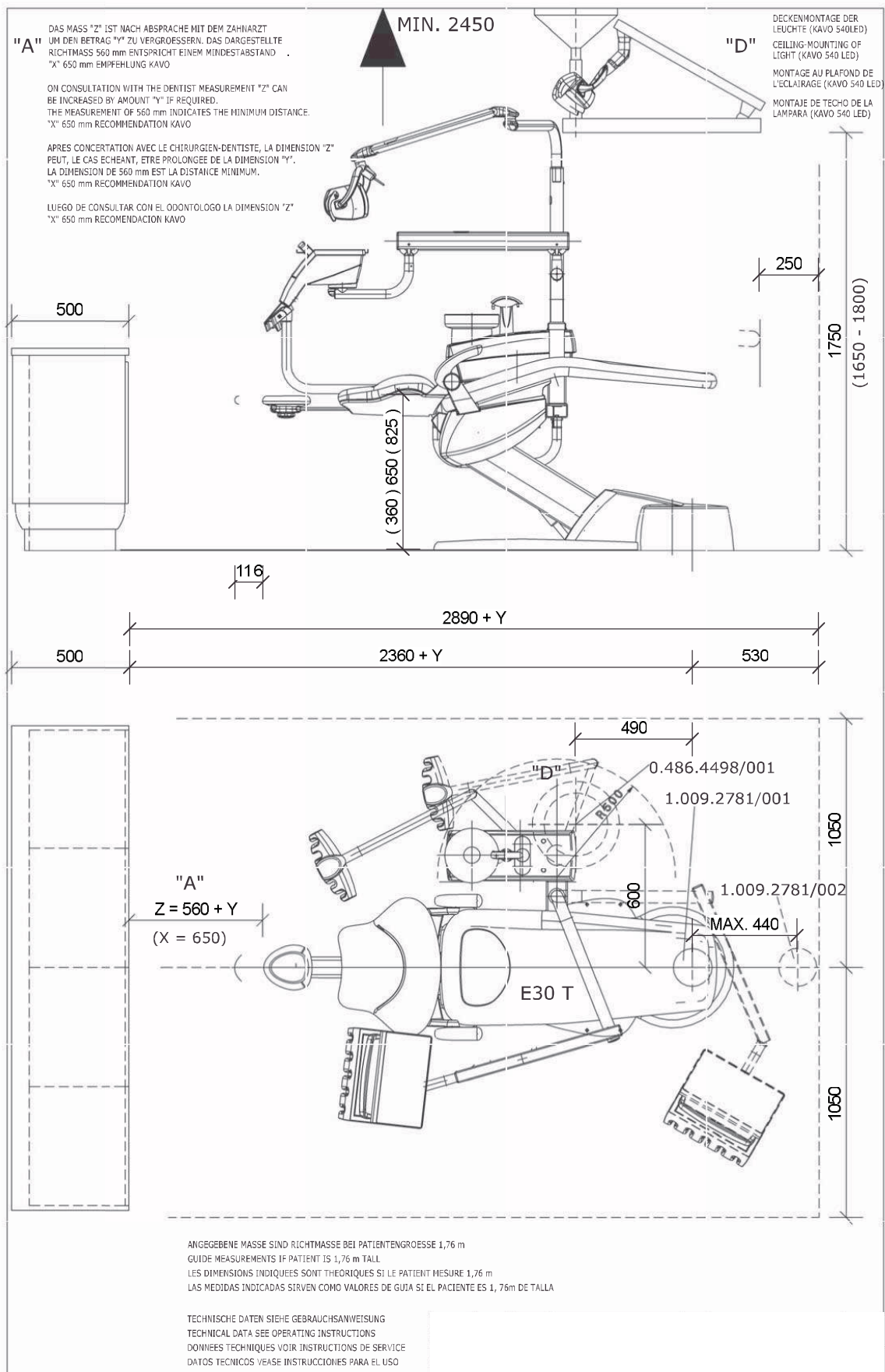

7 Setup plans

